


PolyDetector ALK Magenta Substrate - Chromogen

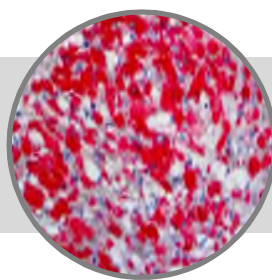

Bio SB
BIOSCIENCE FOR THE WORLD

www.biosb.com

Intended Use

For In Vitro Diagnostic Use.

Summary and Explanation

PolyDetector ALK Magenta Substrate-Chromogen is suitable for use in immunohistochemical, Immunocytochemical and *in situ* hybridization procedures using Alkaline Phosphatase Detection Systems. PolyDetector ALK Magenta Substrate-Chromogen allows for the demonstration of cell antigens or nucleic acids in paraffin-embedded tissues, cryostat sections, cytospins, and cell preparations.

The substrate-chromogen is the final step in the detection portion; it enables the antibody-antigen-enzyme (or nucleic acid-probe-enzyme) complex to be viewed under the light microscope. This occurs when ALK Magenta, in the presence of Alkaline Phosphatase, gets deposited at the site of the target antigen or nucleic acid, producing a magenta color that is partially soluble in organic solvents and therefore care should be taken when mounting with a fast permanent mounting protocol or it can be mounted with permanent medium when using the Bio SB ChromoProtector (BSB 0151 – BSB 0156).

Presentation

PolyDetector ALK Magenta Substrate-Chromogen contains an ALK Magenta chromogen solution, an ALK Magenta Activator solution and an ALK Magenta Buffer-Substrate solution. All the components are buffered, contain stabilizers and a non-azide anti-microbial.

Catalog No.	Buffer - Substrate	Chromogen	Activator
BSB 0077	15 mL	1 mL	1 mL
BSB 0078	50 mL	3 mL	3 mL
BSB 0079	100 mL	6 mL	6 mL
BSB 0080	200 mL	12 mL	12 mL
BSB 0081	1000 mL	50 mL	50 mL

Storage Store at 2-8°C

Stability

This product is stable up to the expiration date on the product label. Do not use after expiration date listed on package label. Temperature fluctuations should be avoided. Store appropriately when not in use, and avoid prolonged exposure to room temperature conditions.

Specimen Preparation

Paraffin sections: This product can be used on formalin-fixed paraffin-embedded (FFPE) tissue sections. Ensure tissue undergoes appropriate fixation for best results. Pre-treatment of tissues with heat-induced epitope retrieval (HIER) is recommended using Bio SB ImmunoDNA Retriever with Citrate (BSB 0020-BSB 0023), ImmunoDNA Retriever with EDTA (BSB 0030-BSB 0033). Tissue should remain hydrated using Bio SB Immuno/DNA Washer solutions (BSB 0029 & BSB 0042).

Frozen sections and cell preparations: This product can be used for labeling acetone-fixed frozen sections and acetone-fixed cell preparations.

Preparation of Working Solution

To prepare a working PolyDetector ALK Magenta Substrate-Chromogen solution, first shake the ALK Magenta Chromogen solution, then add 1 drop of ALK Magenta Chromogen to 1 drop of ALK Magenta Activator and mix. After 10-20, add 1 mL of Substrate-Chromogen solution to the chromogen-activator and mix again. Use this working solution within 5-10 minutes of preparation.

Working ALK Magenta Substrate Chromogen Required	1 mL	2 mL	3 mL
ALK-Magenta Chromogen	1 drop	2 drops	3 drops
ALK-Magenta Activator	1 drop	2 drops	3 drops
ALK-Magenta Buffer-Substrate	1 mL	2 mL	3 mL

Recommended Immunohistochemical Protocol

Step	ImmunoDetector AP	PolyDetector AP	PolyDetector Plus AP
Peroxidase/AP Blocker	5 min.	5 min.	5 min
Primary Antibody	30-60 min.	30-60 min.	30-60 min.
1st Step Detection	10 min.	30-45 min.	15 min.
2nd Step Detection	10 min.	Not Applicable	15 min.
Substrate-Chromogen	5-10 min.	5-10 min.	5-10 min.
Counterstain / Coverslip	Varies	Varies	Varies

Mounting Protocol

a. Alcohol/Xylene Protocol (Fast)

1. After the histological, immunohistochemical or *in situ* hybridization staining procedure is completed, rinse slides in deionized water.
2. Dip the slides in alcohol 30%, 70%, and 100% for 3-5 seconds, then dip for 3-5 seconds in 3 xylenes.
3. Add an organic Permanent Mounting medium such as XyGreen PermaMounter (BSB 0169-0174), PermaMounter (BSB 0094-0097) or similar permanent mounting media.
4. Apply cover slip and air dry before microscopic observation.

b. ChromoProtector Protocol

1. After the histological, immunohistochemical or *in situ* hybridization staining procedure is completed, rinse slides in deionized water. Do not incubate tissue or cell specimens in solvents such as alcohol, toluene, or xylene.
2. Using a coplin jar, or a staining dish with a rack, immerse slides with tissues in ChromoProtector or lay wet slides horizontally and apply sufficient drops of ChromoProtector (BSB 0151 – BSB 0156) to completely cover the tissue. Carefully spread ChromoProtector if needed, but avoid contacting the tissue.
3. Incubate slides for ten minutes at 60 °C to allow ChromoProtector to penetrate tissues.
4. Remove excess ChromoProtector by placing slides vertically over an absorbent material and let excess drain off into absorbent material. Do not rinse slides.
5. Allow slides to COMPLETELY air dry.

NOTE: The ChromoProtector will protect tissue from drying artifacts during the air-drying process.

6. When slides are completely dried, they can be mounted using most standard mounting methods such as aqueous or permanent.

7. Permanent Mounting

- Do not dehydrate slide through alcohol and/or xylene prior to mounting.
- Permanent Mounting medium such as XyGreenPermaMounter (Cat # BSB 0169-0174), PermaMounter (Cat# BSB 0094-0097) or similar permanent mounting media can be added directly to the slide until the tissue or cell specimen is covered.
- If the Permanent Mounting medium does not spread evenly on the dry slide, the slide can be dipped in toluene or xylene for 1 – 2 seconds to aid spreading of the mounting medium.
- Use a minimum amount of mounting medium so that slides dry rapidly.
- Apply coverslip and air dry before microscopic observation.

Product Limitations

Due to inherent variability present in immunohistochemical procedures (including fixation time of tissues, dilution factor of antibody, retrieval method utilized and incubation time), optimal performance should be established through the use of positive and negative controls. Results should be interpreted by a qualified medical professional.


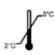





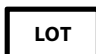
Precautions

1. For professional users only. Results should be interpreted by a qualified medical professional.
2. Ensure proper handling procedures are used with reagent. Minimize microbial contamination of reagents.
3. Always wear personal protective equipment such as laboratory coat, goggles, and gloves when handling reagents.
4. Dispose of unused solution according to local and federal regulations.
5. Do not ingest reagent. If reagent is ingested, seek medical advice immediately.
6. Avoid contact with eyes. If contact occurs, flush with large quantities of water.
7. For additional safety information refer to Safety Data Sheet for this product.
8. For complete recommendations for handling biological specimens please refer to the CDC document, "Guidelines for Safe Work Practices in Human and Animal Medical Diagnostic Laboratories" (see References in this document).

References

1. U.S. Department of Health and Human Services: Centers for Disease Control and Prevention. Guidelines for Safe Work Practices in Human and Animal Medical Diagnostic Laboratories. Supplement / Vol. 61, January 6, 2012.

Symbol Key / Légende des symboles/Erläuterung der Symbole

	EMERGO EUROPE Prinsessegracht 20 2514 AP The Hague The Netherlands		Storage Temperature Limites de température Zulässiger Temperaturbereich		Manufacturer Fabricant Hersteller		Catalog Number Référence du catalogue Bestellnummer
	In Vitro Diagnostic Medical Device Dispositif médical de diagnostic in vitro In-Vitro-Diagnostikum		Read Instructions for Use Consulter les instructions d'utilisation Gebrauchsanweisung beachten		Expiration Date Utiliser jusque Verwendbar bis		Lot Number Code du lot Chargenbezeichnung



69 Santa Felicia Dr., Santa Barbara, CA 93117, USA
Tel. (805) 692-2768 | Tel. (800) 561-1145 | Fax. (805) 692-2769
E-mail: sales@biosb.com | Website: www.biosb.com